




Best™ Radiopaque Gold Marker Strands

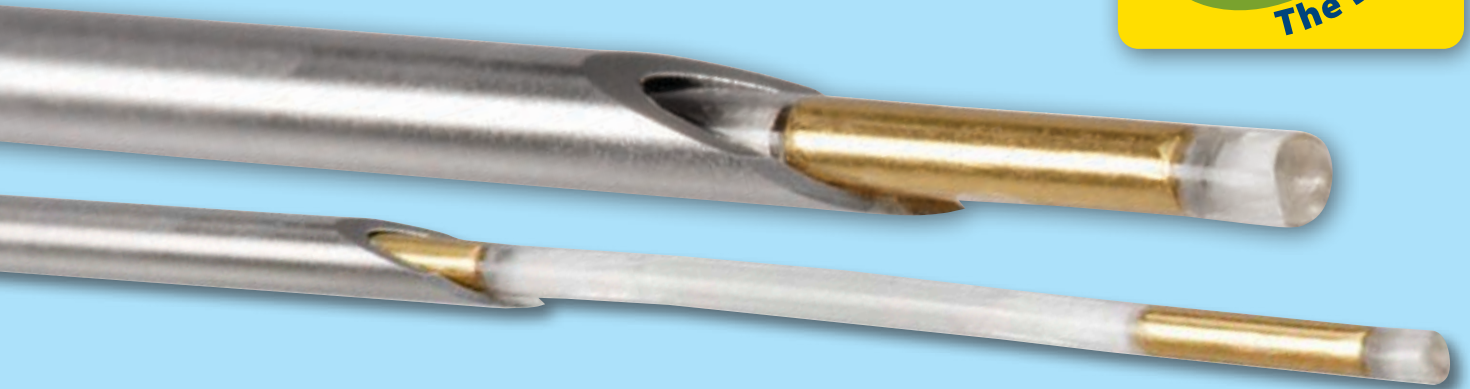

1 cm spaced (center to center)

2 cm spaced (center to center)

0.8 mm x 3.0 mm		0.8 mm x 3.0 mm	
0.8 mm x 5.0 mm		0.8 mm x 5.0 mm	
0.8 mm x 7.0 mm		0.8 mm x 7.0 mm	
1.0 mm x 3.0 mm		1.0 mm x 3.0 mm	
1.0 mm x 5.0 mm		1.0 mm x 5.0 mm	
1.0 mm x 7.0 mm		1.0 mm x 7.0 mm	
1.2 mm x 3.0 mm		1.2 mm x 3.0 mm	
1.2 mm x 5.0 mm		1.2 mm x 5.0 mm	

Stranded single markers

0.8 mm x 5.0 mm	
1.0 mm x 5.0 mm	

Custom strands and loose markers are also available!

All strands pictured have been enlarged from actual size for ease of viewing.

Best Medical International, Inc. © 2019

Best Medical International, Inc. 7643 Fullerton Road, Springfield, VA 22153 USA
 tel: 703 451 2378 800 336 4970 www.besttotalsolutions.com www.teambest.com

AFRICA | ASIA | EUROPE | LATIN AMERICA | MIDDLE EAST | NORTH AMERICA

Best™ Fiducial Marker Kits



Best Medical's Gold Fiducial Markers are available in a variety of diameters and lengths to meet your needs for soft tissue localization such as prostate, liver, pancreas, lung, and breast.

- Best™ Gold Fiducial Markers are available in tissue-absorbable strands
- Standard strand consists of two markers at 2 cm apart center to center
- Special marker strands with custom spacing are also available
- Excellent visualization under X-ray, Fluoro, CT, MRI, Ultrasound and other imaging modalities
- Sold individually, stranded, and in pre-loaded, pre-sterilized needle sets
- Can be used for Brachytherapy, IMRT and IGRT

Best™ Central Rod

These cylinders are to be used with Disposable GYN and Rectal Templates. Cylinders, along with accessories, are now available with the Interstitial Disposable Vaginal Cylinder and a new “Modified” Interstitial Disposable Vaginal Cylinder. The cylinder design, when inserted into the template, positions the Flexi Needles on the surface of the cylinder.



Best™ Central Rod – Modified

The “Modified” cylinder design positions the Flexi Needles 5 mm inward from the surface of the cylinder to reduce dose to the vaginal mucosa. Both designs are available in 15 cm and 25 cm lengths with either a central needle opening or a central tandem opening.




Best™ Multi-Channel Vaginal Applicator




This novel HDR applicator consists of a disposable central insert and an outer sleeve along with a positioning collar. Source channels in the central insert travel fully to the distal end to allow for optimum dose delivery.

Additionally, the outer sleeve provides a 5 mm wall, to reduce dose to the vaginal mucosa, with a 2 mm wall at the distal end of the applicator to allow for greater dose delivery. The applicator is 18 cm in length with diameters ranging from 2.0 cm to 3.5 cm.

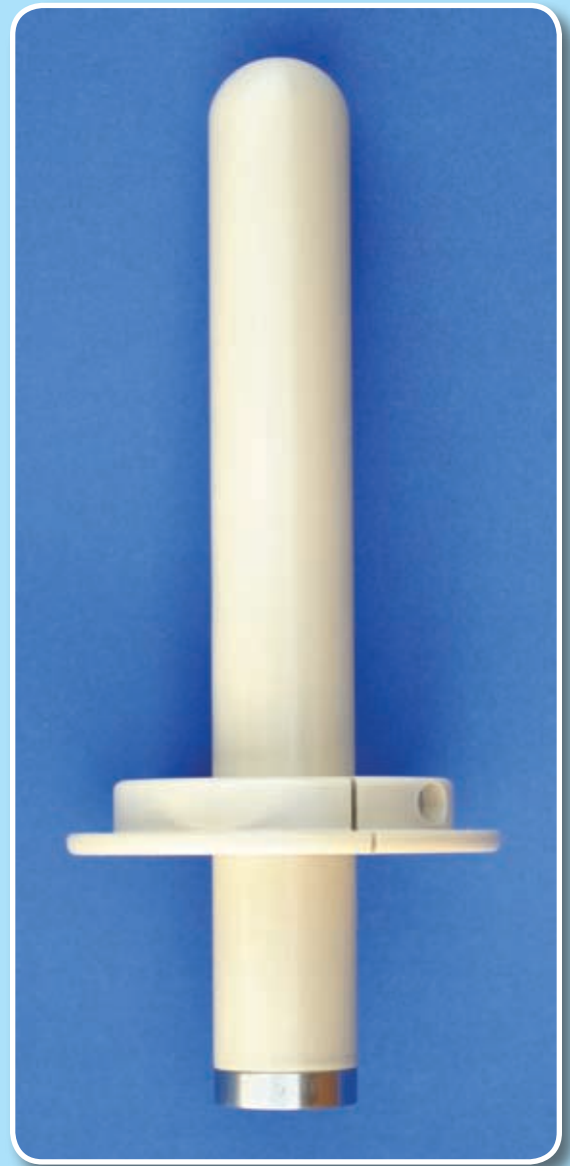


Best™ Brachytherapy Templates


Best Medical offers a wide variety of disposable templates for HDR and LDR applications. Specify hole size when ordering (15, 16 or 17 gauge). These products should be carefully inspected by your staff prior to clinical use and should be used by authorized personnel only.

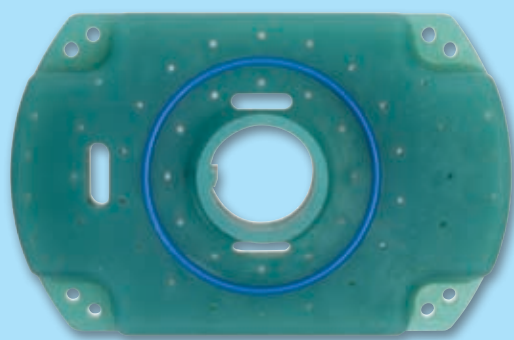


GYN
Template
(Disposable)

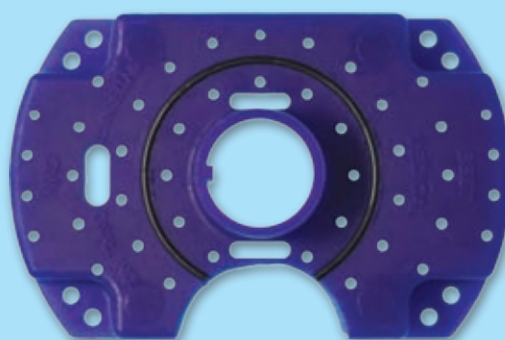


MODIFIED
GYN
Template
(Disposable)

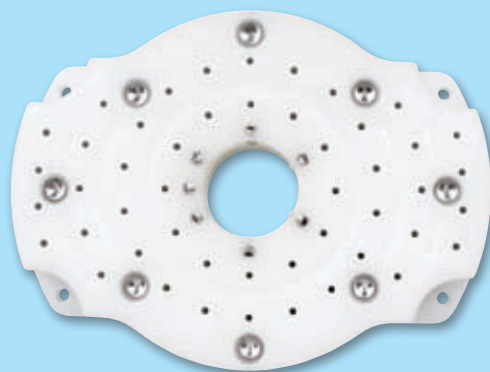


GYN
Template

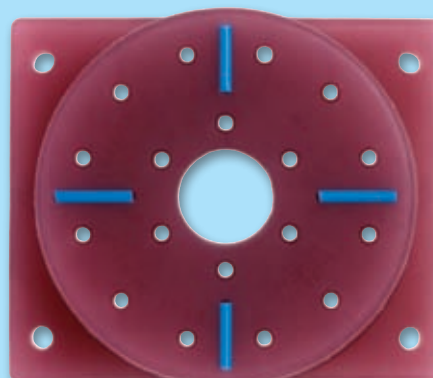

Prostate
Template
(Disposable)

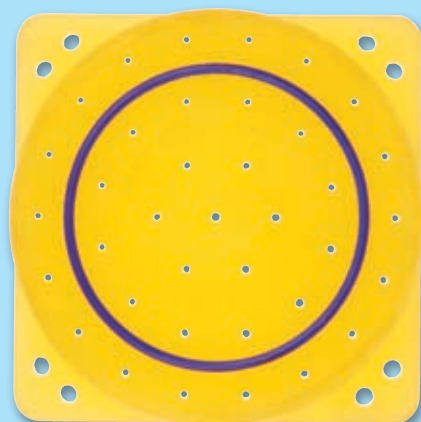


Rectal
Template
(Disposable)

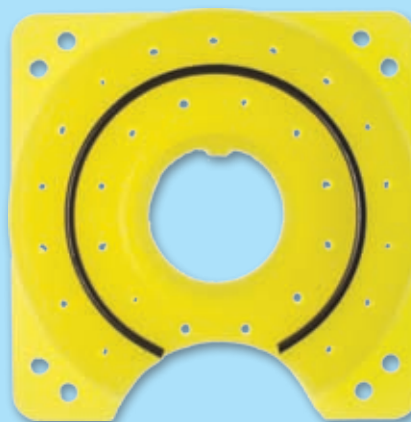

MODIFIED
Rectal
Template
(Disposable)


Best™ X-ray Marker Cables

Best X-ray Marker Cables are indicated for use in HDR brachytherapy procedures. All catheters are available in custom lengths, with or without cable/cable dummy ribbon.

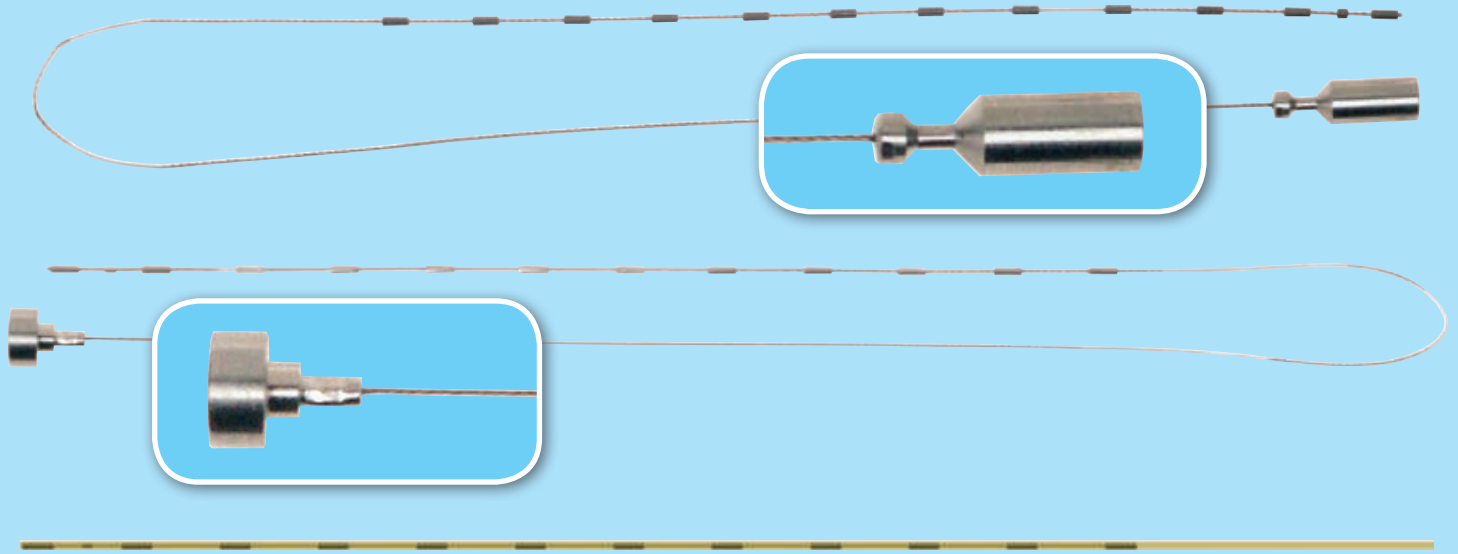


Best™ Prostate Templates

Best Medical International offers a variety of prostate grid templates for all makes and models of steppers, stabilizers and ultrasound systems.

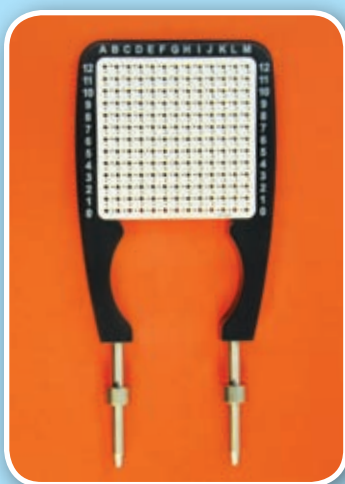

Grid Template (Disposable) with Holder Bracket

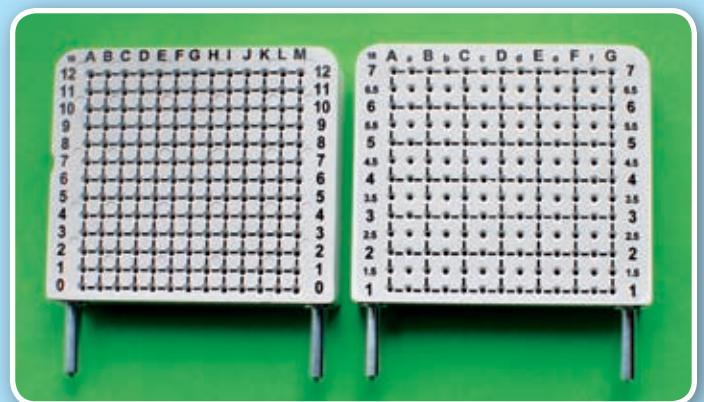

Grid Template for B-K type (Disposable)